


H₂S Scrubbing IsoFlask® with Silver Nitrate

A new option for safe, precise isotopic analyses of sour gases with up to 30% H₂S

This new Iso-brand product delivers the same high-quality, high-precision isotopic results that our customers have come to expect, while addressing concerns associated with the analysis of sour gas:

- ⇒ **The Scrubbing IsoFlask uses a new approach to mitigate the severe isotopic fractionation that occurs with other methods. CO₂ fractionation is completely mitigated with this method.**
- ⇒ **Isotech does not have the capability to work with sour gas samples directly in our facility. The Scrubbing IsoFlask allows the H₂S to be removed by personnel trained to handle H₂S in the field, which allows Isotech to complete the analysis.**

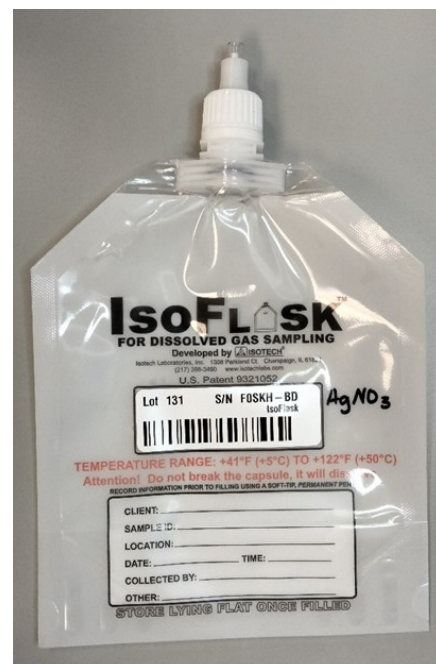

How does it work?

Each Scrubbing IsoFlask comes pre-filled with 20 mL of the scrubbing reagent. The field technician will inject up to 300 cc of sample gas into the IsoFlask and allow the gas to react with the scrubbing solution.

The reaction is fast and visible—the solution will begin to turn black if H₂S is present!

After the reaction has occurred, the field technician will then extract the remaining scrubbed gas from the IsoFlask into an IsoBag, which will be shipped back to Isotech for compositional and isotopic analyses.

Sample Results

Generated by analysis of Isotech reference standard

Component	Composition		δ ¹³ C	
	mol% Expected	mol% Measured	‰ Expected	‰ Measured
Methane	70.0	70.4	-41.00	-40.93
Ethane	10.0	10.1	-28.80	-28.78
Propane	7.07	7.01	-32.37	-32.37
i-Butane	2.54	2.49	-30.98	-30.93
n-Butane	2.55	2.40	-28.14	-28.10
i-Pentane	1.19	1.08	-27.52	-27.45
n-Pentane	1.18	0.985	-30.55	-30.43
CO ₂	5.48	5.37	-40.79	-40.71

Contact Isotech to request a quote!

Call: (217) 398-3490 or E-mail: onlineorders@stratumreservoir.com